

ordering information

Low profile Leo+ baby offers the possibility to treat small arteries down to Ø 1,5 mm.

| | Reference | Vessel Ø (mm) | Stent's length at nominal Ø (mm) | Unconstrained stent | | Compatible delivery catheters |
|---------------|--------------|---------------|----------------------------------|---------------------|-------------|-------------------------------|
| | | | | Ø (mm) | Length (mm) | |
| Leo+ baby 2,0 | LEO. 2,0x12 | 1,5 to 2,5 | 12 | 2,8 | 8 | VASCO+10 |
| | LEO. 2,0x18 | | 18 | | 12 | |
| | LEO. 2,0x25* | | 25 | | 16 | |
| Leo+ baby 2,5 | LEO. 2,5x12 | 2,0 to 3,1 | 12 | 3,2 | 8 | VASCO+10 |
| | LEO. 2,5x18 | | 18 | | 12 | |
| | LEO. 2,5x25 | | 25 | | 16 | |
| | LEO. 2,5x30* | | 29 | | 19 | |
| | LEO. 2,5x35* | | 35 | | 23 | |
| Leo+ baby 3,0 | LEO. 3,0x12 | 2,5 to 3,6 | 12 | 3,8 | 9 | VASCO+10 |
| | LEO. 3,0x18 | | 18 | | 13 | |
| | LEO. 3,0x25 | | 25 | | 16 | |
| | LEO. 3,0x35 | | 35 | | 24 | |
| Leo+ 3,5 | LEO. 3,5x12 | 3,10 to 4,25 | 12 | 4,4 | 9 | VASCO+21 |
| | LEO. 3,5x18 | | 18 | | 12 | |
| | LEO. 3,5x25 | | 25 | | 17 | |
| | LEO. 3,5x30 | | 30 | | 19 | |
| | LEO. 3,5x35 | | 35 | | 24 | |
| | LEO. 3,5x50 | | 50 | | 36 | |
| Leo+ 4,5 | LEO. 4,5x15 | 4,25 to 5,30 | 15 | 5,4 | 12 | VASCO+25 |
| | LEO. 4,5x20 | | 21 | | 15 | |
| | LEO. 4,5x25 | | 26 | | 18 | |
| | LEO. 4,5x30 | | 30 | | 21 | |
| | LEO. 4,5x40 | | 40 | | 28 | |
| | LEO. 4,5x50 | | 51 | | 37 | |
| | LEO. 4,5x75 | | 75 | | 55 | |
| Leo+ 5,5 | LEO. 5,5x25 | 5,30 to 6,50 | 26 | 6,7 | 18 | VASCO+28 |
| | LEO. 5,5x30 | | 31 | | 21 | |
| | LEO. 5,5x35 | | 36 | | 26 | |
| | LEO. 5,5x50 | | 52 | | 35 | |
| | LEO. 5,5x60 | | 60 | | 42 | |
| | LEO. 5,5x75 | | 75 | | 51 | |

*Manufactured on demand.

Reference: 1. Internal data

The self-expandable LEO+/LEO+ Baby stent is designed for the treatment of intracranial aneurysms in association with embolization coils. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT Extrusion SAS. Carefully read the instructions for use before use. First CE marking:2007 (LEO+),2012 (LEO+Baby). VASCO+ is a reinforced micro-catheter intended for injection of diagnostic or therapeutic product, to position pushable coils "SPIRALES" or detachable coils especially the ones of MDS « mechanical detachment system », for the use of the self-expanding stent LEO+ or SILK+. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT Extrusion SAS. Carefully read the instructions for use before use. First CE marking:2004.

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. © 2019 BALT SAS and affiliates, all rights reserved. All representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier's distributors). It cannot be distributed or given to patients. The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. DC015GB (01/2019)

Balt

10, rue de la Croix Vigneron, 95160 Montmorency France

Tel. : +33 (0)1 39 89 46 41

Fax : +33 (0)1 34 17 03 46

www.balt.fr



leo+ & leo+ baby

the authentic braided stent




leo⁺ & leo⁺ baby

Designed for the treatment of intracranial aneurysms in association with embolization coils

aneurysm treatment

More than
40 000
Leos deployed
around the world
since 2003

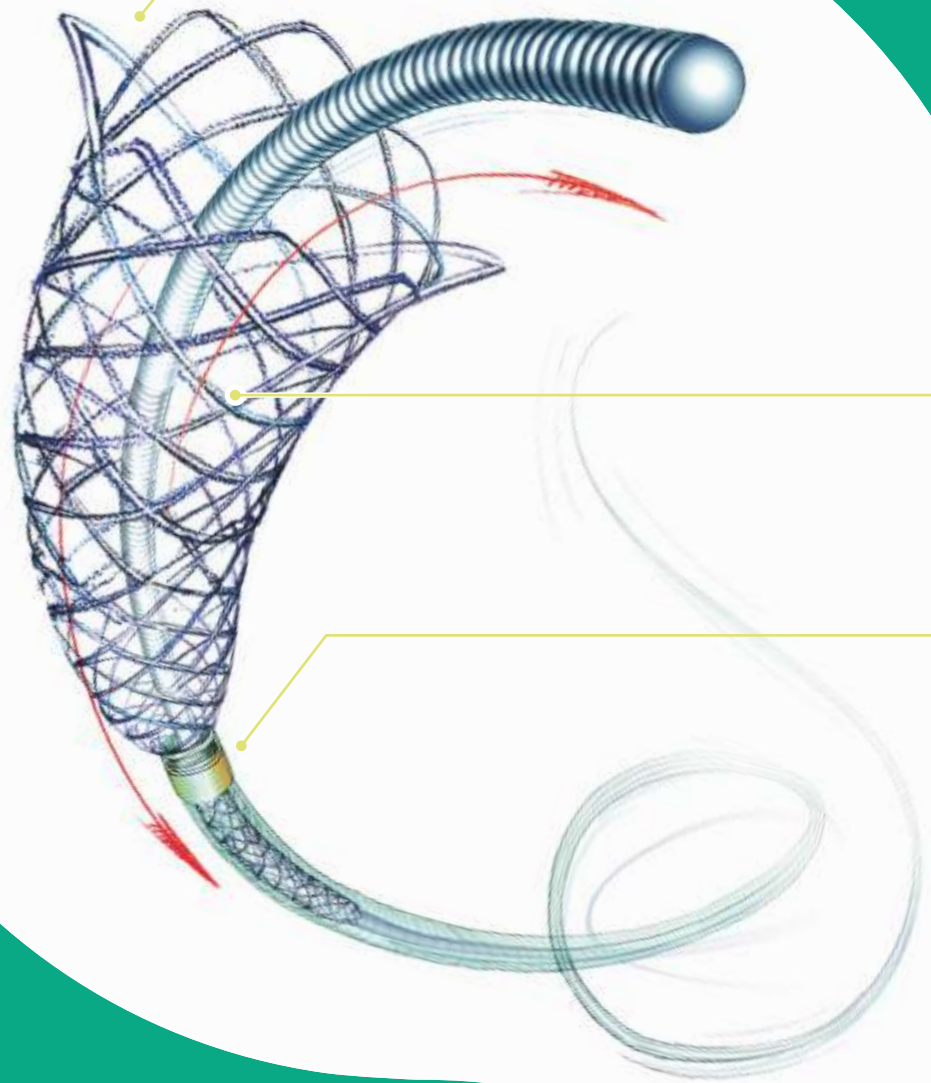

● Customization for
tailor-made solutions

Flexibility
thanks to the nickel titanium materials

Conformability
given by the sliding-cell technology

Smooth wall apposition
as a result of the rounded short flared ends

● Radial force for optimal
coil mass support

Intraluminal support & neck coverage
thanks to the 16 braided wires

● Accurate positioning

Visibility
two helical markers on the
entire body of the stent

Resheathability
up to 90% of deployed length¹

Flow diversion effect

Bouillot et al., Computational fluid
dynamics with stents: quantitative
comparison with particle image
velocimetry for three commercial off
the shelf intracranial stents
J NeuroIntervent Surg 2016;
8:309-315¹

16
braided
wires

resheathable
up to **90%**
of deployed
length¹