

## ordering information

	Reference	Proximal O.D. - Distal O.D.	I.D.	Total length (cm)	Super-supple distal tip length (cm)	Curve
Fargo Mini	FRGMIN4,2F120 FRGMIN4,2F120MP	4,2F to 3,9F	1,02mm (.040")	120	2,5	Straight (D)
	FRGMIN4,2F135 FRGMIN4,2F135MP			135	15	Multipurpose (MP)

	Reference	Proximal O.D. - Distal O.D.	I.D.	Total length (cm)	Super-supple distal tip length (cm)	Curve
Fargo	FRG6F105_8 FRG6F105_8MP	6,0F to 4,9F	1,46mm (.058")	105	8	Straight (D)
	FRG6F115_8 FRG6F115_8MP			115		
	FRG6F125_8 FRG6F125_8MP			125	Multipurpose (MP)	
	FRG6F135_15 FRG6F135_15MP			135		15

	Reference	Proximal O.D. - Distal O.D.	I.D.	Total length (cm)	Super-supple distal tip length (cm)	Curve
Fargo Max	FRGMAX6F95_8 FRGMAX6F95_8MP	6,0F to 6,0F	1,78mm (.070")	95	8	Straight (D)
	FRGMAX6F105_8 FRGMAX6F105_8MP			105		
	FRGMAX6F115_8 FRGMAX6F115_8MP			115	Multipurpose (MP)	
	FRGMAX6F125_8 FRGMAX6F125_8MP			125		

The guiding catheter FARGO are intended to facilitate the introduction of micro-catheters for therapeutic and diagnostic use. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT EXTRUSION SAS. Carefully read the instructions for use before use. Not reimbursed. First CE marking:2009. The self-expandable SILK+ stents are designed for the treatment of intracranial aneurysms and should be used only by clinicians trained in the placement of intracranial stents. Class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT EXTRUSION SAS. Carefully read the instructions for use before use. Not reimbursed. First CE marking:2012. The content of this document, in particular data, information, trademarks and logos is BALT SAS and its affiliates' property. © 2019 BALT SAS and affiliates, all rights reserved. All representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. D:\043GB (08/19)

Balt  
10, rue de la Croix Vigneron, 95160 Montmorency, France  
T l. : +33 (0)1 39 89 46 41  
Fax : +33 (0)1 34 17 03 46  
www.balt.fr

 **balt**  
inspiring innovation

# fargoline

braided guiding  
catheter

access  
treatment



 **balt**  
inspiring innovation

# fargoline

## Braided guiding catheter

Designed to facilitate the introduction of microcatheters for therapeutic or diagnostic use.





wires progressive braiding

**Great stability and support**  
**Back up support, torque & push** provided by a wide braiding in proximal part  
**Kink resistance & navigability** allowed by a tight braiding in distal part

**Navigability through tortuous anatomy**  
**Smooth navigation** thanks to a hydrophilic coated distal part

**Soft & supple** distal part to facilitate navigation in tortuous anatomy

**Gentle access** supported by various distal curves

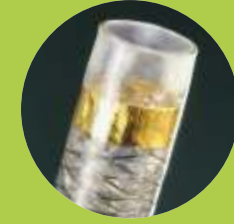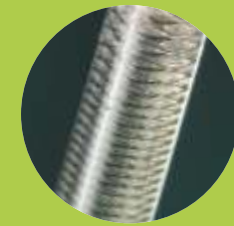

## Large applications

### Fargo mini

#### Triaxial technique

ideal combination between Fargo Max & Fargo Mini providing a high support for lesions in remote anatomies

#### Distal support & stability

to control microcatheter's tip and reduce the catheter's slack

### Fargo

#### More stability & support

ideal for Silk+ cases

### Fargo max

#### Balloon remodeling technique

#### Jailing technique

thanks to its large I.D. of .070"

High flow road map contrast injection

various lengths available

